

Алматы (7273)495-231	Иваново (4932)77-34-06	Магнитогорск (3519)55-03-13	Пермь (342)205-81-47	Тверь (4822)63-31-35
Ангарск (3955)60-70-56	Ижевск (3412)26-03-58	Москва (495)268-04-70	Ростов-на-Дону (863)308-18-15	Тольятти (8482)63-91-07
Архангельск (8182)63-90-72	Иркутск (395)279-98-46	Мурманск (8152)59-64-93	Рязань (4912)46-61-64	Томск (3822)98-41-53
Астрахань (8512)99-46-04	Казань (843)206-01-48	Набережные Челны (8552)20-53-41	Самара (846)206-03-16	Тула (4872)33-79-87
Барнаул (3852)73-04-60	Калининград (4012)72-03-81	Нижний Новгород (831)429-08-12	Саранск (8342)22-96-24	Тюмень (3452)66-21-18
Белгород (4722)40-23-64	Калуга (4842)92-23-67	Новокузнецк (3843)20-46-81	Санкт-Петербург (812)309-46-40	Ульяновск (8422)24-23-59
Благовещенск (4162)22-76-07	Кемерово (3842)65-04-62	Ноябрьск (3496)41-32-12	Саратов (845)249-38-78	Улан-Удэ (3012)59-97-51
Брянск (4832)59-03-52	Киров (8332)68-02-04	Новосибирск (383)227-86-73	Севастополь (8692)22-31-93	Уфа (347)229-48-12
Владивосток (423)249-28-31	Коломна (4966)23-41-49	Омск (3812)21-46-40	Симферополь (3652)67-13-56	Хабаровск (4212)92-98-04
Владикавказ (8672)28-90-48	Кострома (4942)77-07-48	Орел (4862)44-53-42	Смоленск (4812)29-41-54	Чебоксары (8352)28-53-07
Владимир (4922)49-43-18	Краснодар (861)203-40-90	Оренбург (3532)37-68-04	Сочи (862)225-72-31	Челябинск (351)202-03-61
Волгоград (844)278-03-48	Красноярск (391)204-63-61	Пенза (8412)22-31-16	Ставрополь (8652)20-65-13	Череповец (8202)49-02-64
Вологда (8172)26-41-59	Курск (4712)77-13-04	Петрозаводск (8142)55-98-37	Сургут (3462)77-98-35	Чита (3022)38-34-83
Воронеж (473)204-51-73	Курган (3522)50-90-47	Псков (8112)59-10-37	Сыктывкар (8212)25-95-17	Якутск (4112)23-90-97
Екатеринбург (343)384-55-89	Липецк (4742)52-20-81		Тамбов (4752)50-40-97	Ярославль (4852)69-52-93

Россия +7(495)268-04-70

Казахстан +7(7172)727-132

Киргизия +996(312)96-26-47

www.ami.nt-rt.ru | | aix@nt-rt.ru

Технические характеристики на гидропилингвые системы Aqua Sonic компании AMI INC.

Aqua Sonic

The Newest wonderful peeling for Skin Resurfacing

Non-Invasive

-

Painless

-

Multifunctional

L. J



Technical Specifications

Ultrasonic	1.2 MHz
Motor	upto 600mmHg (±10%)
Power consumption	80Watt
Input Voltage	220VAC / 50 / 60Hz
Dimensions	350(D) X 350(L) X 1000 (H) mm
Weight	30Kg
Aqua Hand piece	Face(Regular Size) Nose(Small Size)
Aqua Waste tank	1200cc

Skin Refreshing!

Make your skin cleaner and whiter.

Hydro Peeling System

Aqua Sonic

Remove horny substances and supply moisture at the same time.

Aqua Sonic is designed to promote exfoliation to remove aged dead skin cells and sebum the same time, injects highly concentrated nutritious serum to make skin and scalp moist, smooth, and lustrous.

Safe, no downtime, no side effect.

Because the treatment does not cause any damage to skin and patients can lead a normal life just after treatment. Also there is no side effect such as injury, pain, burns, rash, etc.

Multiple skincare effects.

Treatment removes foreign materials such as cosmetics and sebum. Therefore, it is expected to achieve whitening and pore-tightening effects.



Aqua Sonicis...

Not a decortication old style machine. It is a new-concept natural spa skincare system, which makes your skin exfoliation with no pain. Also, a range of treatment is variety such as cleansing of skin and removal of obsolescent horny substances. If you use Laser Doctor as a helper, it makes your skin nutritious, calm and Epidermal Growth Factor (EGF). After all the treatment, you will feel moist and refreshing. Ultra Sound can release pain, tighten pores and treat for fat.

Name and function of parts and components

Aqua Sonic is designed to promote exfoliation to remove aged dead skin cells, sebum and hydrate your skin. Laser Doctor as a helper makes your skin nutritious, calm and Epidermal Growth Factor (EGF).

Because treatment does not cause any damage to skin and scalp, patients can lead a normal life right after treatment. There are no side effects such as injuries, pain, burns, rashes, etc.

Treatment removes foreign materials (such as cosmetics and sebum) so users can expect to achieve whitening and pore-tightening effects.



Aqua Sonic adopts a vacuum extrusion method. The hand piece should come into close contact with the skin. You can control pressure for both fine bubbles and serum liquid to perform treatment. Red areas on the skin that may be created during operation will disappear normally in 20-30 minutes and the skin restores healthy and lustrous condition.



Wipe wastes 2-3 times off the surface of the face with a cleansing cream

Diagnosis of skin and Selection of serum

Identify whether the skin is normal, oily, sensitive, or of acne type and make selection of a serum suitable for the skin type. [Aqua PE: For all skin type]

Treatment [removal of keratin and sebum]


First spread the selected fine bubble and serum over the face and then divide the face into 3-5 areas to start treatment for about 2 minutes.

Washing

Wash and wipe out the face cleanly with a spongy after treatment.

Treatment (supply of moisture and nutrients) and finish

Perform treatment while letting serum penetrate into the skin; after treatment, give a finishing touch with a skin balm or a nourishing cream.

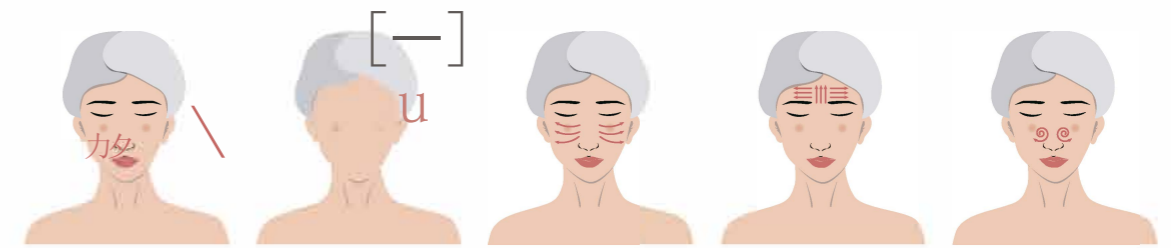

In the order of Neck Chin Cheek Forehead Nose

- 1 Neck after turning face a little sideway, make lymph nodes relax by massaging the area from top to bottom and in the opposite way with the sternocleidomastoid region in the center.
- 2 Nose perform treatments repeatedly because there is lots of sebum in this region
- 3 Perform treatment by traveling all across the face from bottom to top along the texture of muscles
- 4 Treatment is recommended to be made for about 2 minutes respectively in the order of the neck chin cheek forehead nose

Aqua Sonic Solution

For removal of keratin as well as for extrusion



The special tip of the Aqua Sonic consists of two suction holes and one serum discharge hole. Its specially-designed spiral tip is intended to hold a product long enough to soften dead cells of skin, provide nutrients, and absorb them in a vacuum state to maximize its peeling and extrusion effects.



AQUA PE [For all skin Type- For the general skin, 30ml]

Since glucosamine and lactic acids promote the peeling of the skin to remove impure substances softly and plant extracts soothe and moisten the skin smoothly, the product enables a scientific treatment that does not produce any injuries or rashes while achieve a bright and healthy-looking skin regeneration.

Dr. Serum [Nutrition and supply of moisture for all type of skin, 20ml]

The deep penetration of vitamin C, E, Aloe Vera, and green tea extract will help skin to restore and regenerate.

RF Cream [Lifting and Slimming]

This cream is formulated to use with High Frequency equipment. It is well rubbed until the treatment finish. It helps grow and maintain healthy skin with the botanical extract and anti-oxidant ingredients such as Allantoin, Green Tea Extract, Portulaca Extract, Belamcanda Chinese Water Extract, Tocopheryl Acetate (anti-oxidation ingredients) and Beta Glucan.

Алматы (7273)495-231
Ангарск (3955)60-70-56
Архангельск (8182)63-90-72
Астрахань (8512)99-46-04
Барнаул (3852)73-04-60
Белгород (4722)40-23-64
Благовещенск (4162)22-76-07
Брянск (4832)59-03-52
Владивосток (423)249-28-31
Владикавказ (8672)28-90-48
Владимир (4922)49-43-18
Волгоград (844)278-03-48
Вологда (8172)26-41-59
Воронеж (473)204-51-73
Екатеринбург (343)384-55-89

Иваново (4932)77-34-06
Ижевск (3412)26-03-58
Иркутск (395)279-98-46
Казань (843)206-01-48
Калининград (4012)72-03-81
Калуга (4842)92-23-67
Кемерово (3842)65-04-62
Киров (8332)68-02-04
Коломна (4966)23-41-49
Кострома (4942)77-07-48
Краснодар (861)203-40-90
Красноярск (391)204-63-61
Курск (4712)77-13-04
Курган (3522)50-90-47
Липецк (4742)52-20-81

Магнитогорск (3519)55-03-13
Москва (495)268-04-70
Мурманск (8152)59-64-93
Набережные Челны (8552)20-53-41
Нижний Новгород (831)429-08-12
Новокузнецк (3843)20-46-81
Ноябрьск (3496)41-32-12
Новосибирск (383)227-86-73
Омск (3812)21-46-40
Орел (4862)44-53-42
Оренбург (3532)37-68-04
Пенза (8412)22-31-16
Петрозаводск (8142)55-98-37
Псков (8112)59-10-37

Пермь (342)205-81-47
Ростов-на-Дону (863)308-18-15
Рязань (4912)46-61-64
Самара (846)206-03-16
Саранск (8342)22-96-24
Санкт-Петербург (812)309-46-40
Саратов (845)249-38-78
Севастополь (8692)22-31-93
Симферополь (3652)67-13-56
Смоленск (4812)29-41-54
Сочи (862)225-72-31
Ставрополь (8652)20-65-13
Сургут (3462)77-98-35
Сыктывкар (8212)25-95-17
Тамбов (4752)50-40-97

Тверь (4822)63-31-35
Тольятти (8482)63-91-07
Томск (3822)98-41-53
Тула (4872)33-79-87
Тюмень (3452)66-21-18
Ульяновск (8422)24-23-59
Улан-Удэ (3012)59-97-51
Уфа (347)229-48-12
Хабаровск (4212)92-98-04
Чебоксары (8352)28-53-07
Челябинск (351)202-03-61
Череповец (8202)49-02-64
Чита (3022)38-34-83
Якутск (4112)23-90-97
Ярославль (4852)69-52-93

Россия +7(495)268-04-70

Казахстан +7(7172)727-132

Киргизия +996(312)96-26-47

www.ami.nt-rt.ru | | aix@nt-rt.ru